

ABOUT US - MISSION

HAROL is the first Italian company that over the past 30 years, thanks to close collaboration with doctors and nurses, has developed various models of technologically advanced interfaces for the patient.

These solutions fully reflect and satisfy the evolution of the operators' needs: **high quality, high performance and maximum reliability combined with compactness and portability, without ever losing sight of patient comfort.**

HAROL

Innovative solutions

We have been on the market for several years with the COMPACT "family".

Compact, practical and adaptable devices for any emergency situation for intra and extra-hospital use.

They allow CPAP therapy even with a single oxygen source, by delivering high flows.


GENERAL FEATURES

	COMPACT I - HAR	COMPACT II - HAR	COMPACT III - HAR
TYPE OF MASK AND HEADGEAR	Mask with inflatable cushion and disposable headgear	Mask with inflatable cushion and disposable headgear	Mask with inflatable cushion and disposable headgear
AVAILABLE SIZES	5 sizes XS - S - M - L - XL and color code	5 sizes XS - S - M - L - XL and color code	4 sizes S - M - L - XL
OXYGEN SUPPLY TUBING FOR STANDARD OPERATIONS	Single O ₂ line (port for second line)	Single O ₂ line (port for second line)	Double O ₂ line
FLOW-METERS FOR STANDARD OPERATIONS	A common 15 LPM (3.5 bar) flow-meter	A common 15 LPM (3.5 bar) flow-meter	Double 30+30 LPM (3.5 bar) flow-meter
FiO ₂ % CONCENTRATION	From 40 to 60 O ₂ with single line. Integrated manometer for measuring PEEP	From 50 to 60 O ₂ with single line. Integrated manometer for measuring PEEP	Adjustable FiO ₂ from 30 to 100. Integrated manometer for measuring PEEP
PEEP VALVE	Adjustable from 5 to 15 cmH ₂ O	Virtual from 3 to 15 cmH ₂ O	Adjustable from 0 to 20 cmH ₂ O
NEBULIZER FOR MEDICINES AND MDI PORT	YES	YES	YES
HEPA FILTER	YES	YES	YES

THE BEST SOLUTIONS FOR EMERGENCY CPAP

COMPACT I - HAR

- » **COMPLETE SYSTEM**
with one-way facemask.
- » **COMPACT VENTURI DEVICE**
works with a single oxygen source.
- » **INTEGRATED PEEP VALVE**
adjustable from 5 to 15 cmH₂O.
- » **MANOMETER**
for measuring PEEP, with integrated anti-suffocation valve.
- » **ADJUSTABLE FiO₂**
from 40 to 60.
- » **POSSIBILITY**
of controlling the HAR-NEB nebulizer system (optional accessory) without the use of an additional oxygen source.
- » **MDI PORT**
for the administration of medicines.
- » **ANTIBACTERIAL FILTER**
for the protection of the patient and health professionals.

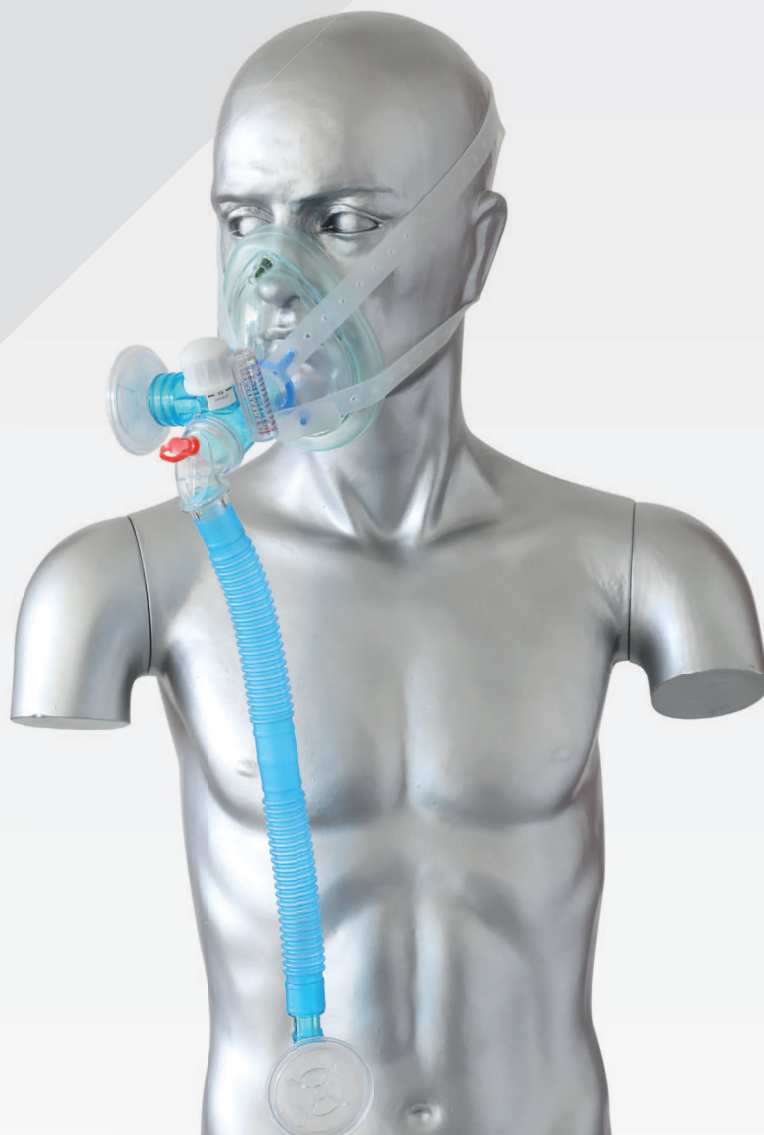

TECHNICAL FEATURES

Set. CPAP	5 cmH ₂ O	6 cmH ₂ O	7,5 cmH ₂ O	10 cmH ₂ O	12,5 cmH ₂ O	15 cmH ₂ O
Supply Oxygen	5 LPM	6 LPM	7 LPM	8 LPM	9 LPM	10 LPM
FiO ₂ %	39 %	43 %	50 %	55 %	57 %	59 %

AVAILABLE PRODUCTS

ITEM NUMBER	SIZE	PACKAGING
CPAP10208	PEDIATRIC	Single packed
CPAP10209	SMALL FOR ADULTS	Single packed
CPAP10210	MEDIUM FOR ADULTS	Single packed
CPAP10211	LARGE FOR ADULTS	Single packed
CPAP10212	EXTRA LARGE FOR ADULTS	Single packed
CPAP10700	VENTURI DEVICE	Single packed